


medivon.

User manual

Professional Air
Massage Mat

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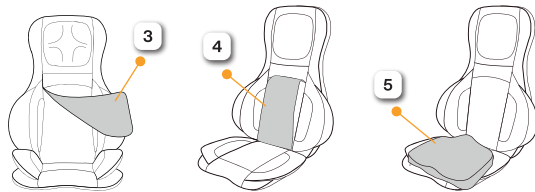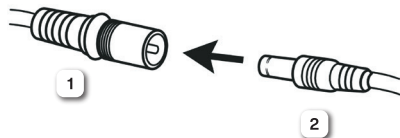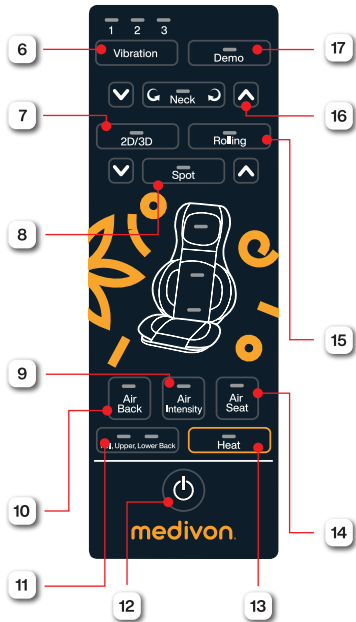1. Configuration and usage

- Attach the mat to the chair with the elastic straps
- Connect Mat connector (1) with the power supply connector (2)
- Plug the home power adapter into an electrical outlet
- Turn on the device with the remote control (instructions later)
- Place the RC in the bag on the side of the massage mat

Note: The removable flap (3) or cover for the neck massage section allows you to adjust the intensity between a more gentle or more intense massage. If the massage is too strong, place a blanket or towel between you and the mat (4)

If you want to massage your neck, shoulder, lower back, but the shiatsu massage heads cannot reach these areas, you can adjust the height of the massage by placing a pillow on the sitting area (5).

To avoid overloading the muscles, we recommend that you do not exceed 15 minutes for one massage session. Do not massage yourself more than 3 times during the day. (in the full range of device operation - 15 minutes)

2. Care and cleaning

- Do not machine wash or submerge the product or adapter in water.
- Always turn off the power before cleaning and remove the adapter from the power supply with clean, dry hands.
- Wipe the device with a slightly damp cloth and dry it thoroughly.
- Do not wipe the device with preparations containing alcohol or harsh chemicals such as thinner or bleach.
- The neck cover is removable and washable.
- Do not stretch the elastic straps excessively on the back.

- Store away from heat, fire or direct sunlight.
- Do not put excessive pressure on the device and the massage heads. Applying a direct force in excess of 49 kg may damage the massager and void the warranty.
- Store the massager in a dry place inaccessible to children.
- Avoid contact with sharp edges or sharp objects that could cut or puncture the material of the device.

3. Remote controller description

- (6) - press to turn on the massager vibration, you can switch between 3 levels of intensity.
- (7) - press to change the 2D / 3D massage mode. In the 3D option, it additionally takes into account the exact area of the spine.
- (8) - press to keep the shiatsu function in one area. Press the arrow buttons to move the shiatsu heads up or down.
- (9) - press to set the intensity of air massage on low, medium or high level.
- (10) - press to enable or disable rear air massage.
- (11) - press to designate the area with shiatsu massage.
- (12) - press to turn the device on or off.
- (13) - press to turn the heating on or off. The heating cannot be turned on if the back shiatsu function is not turned on.
- (14) - press to enable or disable seat air massage.
- (15) - press to turn on the roll massage function. Repeated pressing changes modes.
- (16) - press to start shiatsu massage on the neck. The arrows change the position of the massage heads in the mat.
- (17) - a quick demonstration of all the functions of the massage mat

4. Safety instructions

- The device should never be left unattended when connected to the power supply. Before assembling or disassembling, disconnect the device from the power supply.
- Do not use the device under a cover where the ventilation opening may be blocked. The device may overheat.
- This device requires close supervision by adults when it is used by children or people with disabilities.
- Do not use with other accessories not recommended by the manufacturer.
- Never use the appliance if it has a damaged cord or plug, if it does not work properly, if it has been dropped, damaged or dropped into water. In such cases, contact the service immediately.
- Do not carry the device with the power cord or use it as a handle.
- Never use the device with blocked air openings.
- Never insert any objects into any openings.
- Do not use outdoors.
- When using, do not wear jewelry or loose clothing.
- Longer hair must be tied back before use so that it does not come into

contact with the device.

- Do not use the warm-up function in places where you have skin disorders, open wounds, swelling or inflammation.
- If the plugs or cables are damaged, they must be replaced by the customer service
- In case of failure, disconnect the device from the mains immediately.
- Do not use in places if you have skin abnormalities, open wounds, swelling or inflammation.
- Improper use or incorrect use excludes all liability for damages.
- Do not use this device while driving or sleeping.
- We have equipped the device with overheating protection for safe use.
- Do not use the device during pregnancy.
- Do not use if you have pacemakers or are under treatment for diseases or other ailments.

The crossed out garbage can symbol means that in the European Union, after using the product, it should be disposed of at a separate, dedicated point. Do not dispose of these products with unsorted municipal waste.

MEDIVON sp. z o. o. Hereby declares that this massage device complies with the European Directive RE 2014/35/EU The full text of the Declaration of Conformity is available at the following internet address:
www.medivon.pl/declarations

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