
medivon.

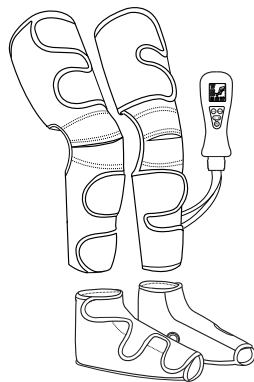
## **Air massager**

Pressure therapy

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## 1. Panel Sterowania

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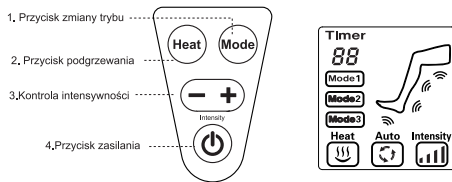
## 1. Control Panel
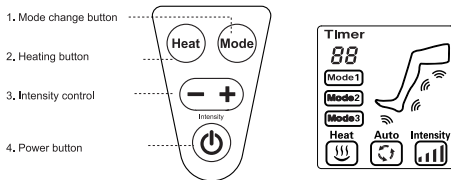
## 2. Operation of the device

 Press this button to turn on the machine when it is in standby. The device will start up with the LCD screen and backlight. The default state is automatic air compression mode, intensity level 2 and heating level 2 - the higher level. A countdown will begin (default timer: 20 minutes). The LCD display shows the current operating mode, the intensity and the remaining operating time. If the button is pressed when the device is turned on or after 20 minutes of device operation, all massage operations will be stopped and the device will be turned off and go into standby mode. When the device is turned on, the backlight will be automatically turned off 15 seconds after the last correct button press. When the backlight is off, press any button to turn it on and all function buttons will be active.

Press this button when the device is working to switch the pressure massage modes. Each time the button is pressed, the massage mode switches respectively once. The massage intensity is automatically adjusted to level 2. The LCD screen shows the massage mode and intensity. The switching sequence is: Auto> Mode 1> Mode 2> Mode 3> Auto. Mode 1 massages only the foot area, Mode 2 massages only the calf area, Mode 3 massages only the thigh area, and Auto Mode massages the entire leg.

Press this button when the appliance is running to change the status of the heating function. Each time the button is pressed, the heating status changes once. The  icon is displayed when the heating function is on. The  icon flashes at heat level 1 and lights up at heat level 2.

Press this button to increase or decrease the massage intensity. Each time the button is pressed, the air compression intensity changes to a higher or lower level, up to level 4. The corresponding intensity icon will be shown on the LCD display.

This icon indicates the area of the body that is currently being massaged

### 3. Dangers

1. Do not use spare parts that are not approved by the manufacturer.
2. Unprofessional service should not disassemble and replace parts of the leg massager without permission.
3. If the power cord is damaged, it must be replaced by appropriate professional personnel to avoid a hazard.
4. Do not use or store the product in a place with high humidity or high temperature.
5. After completing the massage, disconnect the power cord.
6. Do not use the product for more than 15 minutes at a time.
7. If any side effects occur, please stop using the massager immediately.
8. If the temperature is too high or the product starts to smell bad or smells burnt, stop using it immediately and disconnect the power.
9. Do not carry the product by pulling on the power cord to prevent the power cord from loosening.
10. This product is not intended for people who have no physical, sensory or mental abilities or have no experience and knowledge in this field (including children).
11. People insensitive to heat should use the heating function carefully - preferably under the supervision of a caregiver.

The crossed out wheellie bin symbol means that in the European Union, after using the product, it should be disposed of at a separate, dedicated point. Do not dispose of these products with unsorted municipal waste.

### 4. Warnings

If you have any of the following conditions, consult your healthcare practitioner before use:

- Pregnant women, children or the elderly
- People who suffer from skin infections or have reddened and burned skin
- who use a pacemaker or similar implanted medical device;
- Patients with heart disease or malignant tumors, etc.
- People who are treated or who experience discomfort
- People suffering from osteoporosis
- People with high blood pressure or abnormal blood pressure

### 5. Cleaning and maintenance

Turn off the product and unplug the power cord before cleaning. Wipe the surface of the product with a dry, soft cloth. Do not clean the product and accessories with harsh detergents such as benzene, thinners, and store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or humid environment, otherwise it may cause color fading or product damage.

MEDIVON sp. z o. o. Hereby declares that the massage device complies with the European Directive RE 2014/35 / EU. The full text of the Declaration of Conformity is available at the following internet address: [www.medivon.pl/deklaracje](http://www.medivon.pl/deklaracje)

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